# *pearl*

HAIR REMOVER

Model / Modelo PH01

# WELCOME TO PEARL!

You have just purchased a new at-home professional hair-removal system that is completely safe and painless, and that will become your best ally for eliminating unwanted hair.

Pearl is a new experience in hair removal. It is different from the traditional methods you are familiar with.

With Pearl's hair-removal heat technology, heat is transmitted to the hair follicle, which weakens it and reduces hair density after a few months of continuous use.

Until you get used to using the product, it may take 30 to 45 minutes to remove the hair on your legs to begin with, but you'll master the gliding motion and become an expert within just two or three sessions.

Very soon you'll be able to reduce how often you use Pearl because the hair growth will slow down with every use.



**Pearl: long-lasting, painless hair removal**


Before using Pearl for the first time, please read the user manual from start to finish.

# CONTENTS

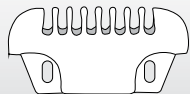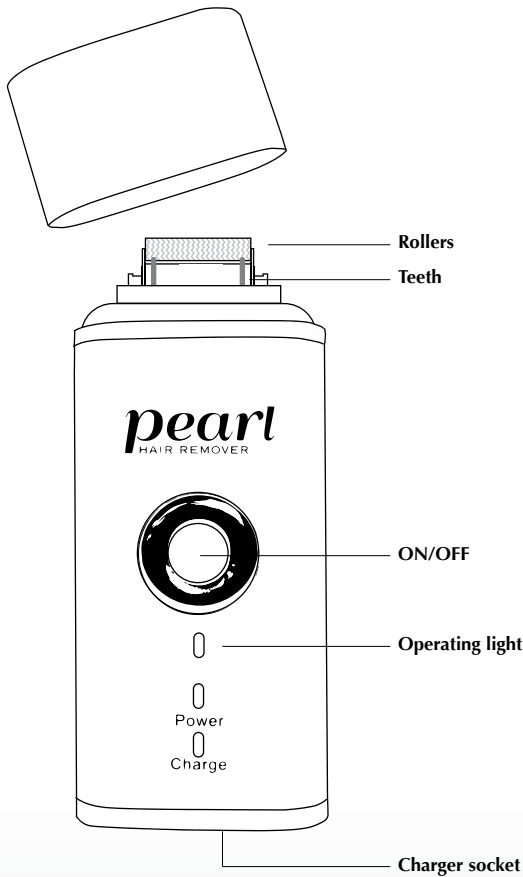
Large head

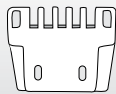Small head


Main unit

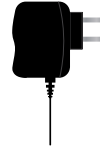
1 charger

2 heads

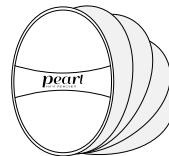Set of buffer pads

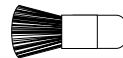1 cleaning brush

1 travel bag

## CHARGING

Charge the Pearl device for five hours before using it for the first time. The green "charge" light will appear to indicate that the battery is charging. The green light will go off to indicate that the device is fully charged. When the device has finished charging, insert the head to begin using it.

**Important:** In order to check the charge status of the battery, the head must be inserted. For safety reasons, it is not possible to switch Pearl on if the head is not inserted.

The charge will last approximately 45 minutes and the battery will take approximately 90 minutes to recharge after the first use. Pearl is equipped with a safety feature that causes the device to switch off automatically if it has not been in use for 15 minutes.

**Important:** Ensure that your skin is clean, dry and free from lotions and oils before using Pearl.

## REMOVE THE COVER AND INSERT THE HEAD

Choose the appropriate head: large for areas such as the arms, legs and back, or small for more delicate areas such as the face, armpits and bikini line.

Line up the pins of the device with the holes in the head and gently push it into place.

## SWITCHING IT ON


Hold the ON/OFF button down for one second to switch the device on.

## GLIDING

- Hold the head of the device against the skin to form a 90° angle with the area to be treated. Start gliding and hold it with the roller flat against the skin.


- Continue gliding and applying pressure on the skin using a smooth, constant, upward movement against the direction of hair growth.

- A red light will come on and stay on as long as you maintain the correct angle and speed during use.

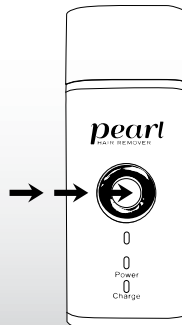- **Important:** Keep the red light steady. The red light will flicker if you do not maintain the correct gliding movement. Adjust the speed and/or angle accordingly until the red light is steady.

- The smell of burning hair and even the appearance of a small amount of smoke means that Pearl is working correctly by crystallizing the hair on contact. A scented candle or air freshener may be used to help mask the smell.
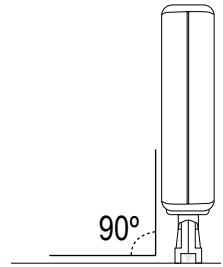
### USING THE BUFFER PAD

To achieve a smooth finish, firmly massage the treated skin back and forth or in circular movements to remove any traces of unwanted hair that have crystallized on the skin.

Apply moisturizing cream to the area after each treatment to achieve a smooth, even finish.

### CARE AND MAINTENANCE



Switch off the device and wait for five minutes before cleaning. Always disconnect the device before starting the cleaning process.

After each use, use the cleaning brush to remove any hairs left behind in the head, taking care not to damage the head or the wire.

### REPLACEMENT OF THE HEADS

After repeated, long-term use, the wire inside the head will wear out and may sometimes break. This is completely normal and to be expected.

Do not attempt to repair it or reuse a head that has broken or become damaged in any way. Discard it and replace it with a new head each time this happens.
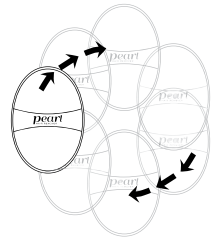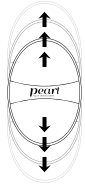
## ADVICE FOR USE

- For best results, the hair to be removed should be between 1 and 3 mm long. If the hair is too short, it will be difficult to remove, since Pearl is designed to avoid direct contact with the skin and therefore prevent risk of burns. If the hair has become ingrown, the skin should be dampened with hot water and then exfoliated. This can be done using exfoliating cream or a loofah mitt. Use circular movements, starting at the ankles and moving up towards the thighs. This ensures that the pores open and makes hair removal more effective.

- Ensure that the head is always clean and dry.

- Choose the appropriate head: large for areas such as the arms, legs and back, or small for areas such as the face, armpits and bikini line.

- In concave areas or areas with folded skin, stretch the skin and make sure it is as smooth as possible to prevent accidents.

- Follow the natural curve of the body when treating rounded areas such as the knees, elbows and ankles.

- **IMPORTANT:** Use the buffer pad to achieve a smooth, even finish. A moisturizing cream or lotion should also be used after treatment. These steps are essential to guarantee satisfactory results.

- Ensure that the heads are always kept clean.

- For the first four to six weeks, use Pearl as many times a week as necessary.

- After the first four to six weeks and as time passes you will notice that you can reduce how often you use Pearl.

## ADVICE FOR TREATING FACIAL HAIR

Read these important recommendations for a more effective treatment.

### TEST

Before you begin, read the instructions and treat a small, discreet patch on your face to see how your skin reacts to the facial hair treatment. If you experience any discomfort, **STOP**.

### UPPER LIP AREA

Start from the corner of the mouth and work upwards, taking care not to touch the lips. Hold the skin tightly and/or purse your lips inwards to create the smoothest surface possible for the gliding action of Pearl. **DO NOT** allow the lips to come into contact with the wire inside the head. **DO NOT** go over one area excessively. If you experience any discomfort, **STOP**.

### UPPER AND JAW LINE


Working from just below the chin or jaw line, glide Pearl around the curve. If necessary, pull the skin tight and/or purse your lips inwards to give your chin and jaw a straighter, smoother line. If you experience any discomfort, **STOP**.

### SIDEBURNS

Start at the base of the chin and glide upwards, following the curve of your face. If necessary, pull the skin tight and straight for a smooth, flat working surface. If you experience any discomfort, **STOP**.

### IMPORTANT NOTE:

The area around the eyes is especially sensitive. Do not treat this area.

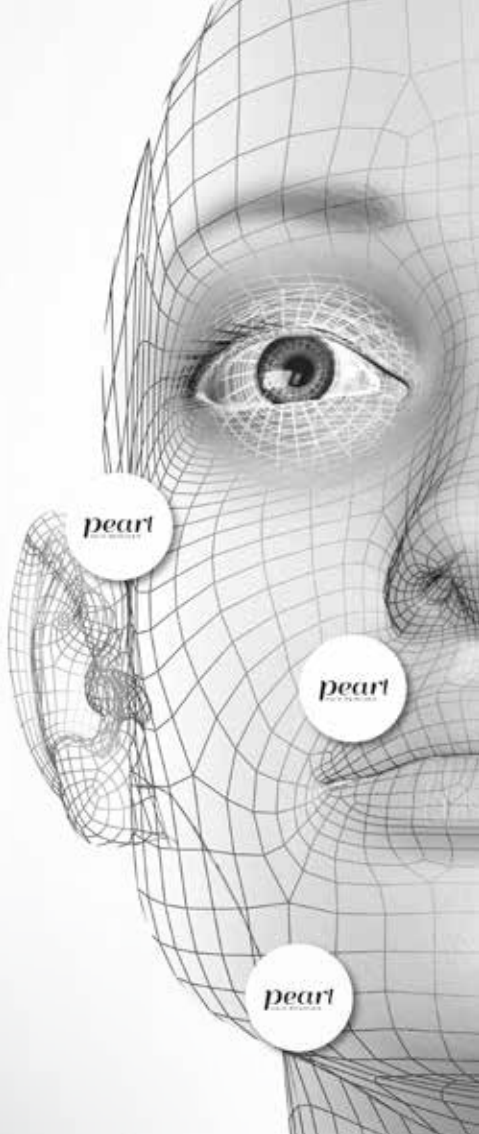**SAFETY INSTRUCTIONS**  
**READ ALL INSTRUCTIONS BEFORE USE.**  
**ALL WARNINGS AND PRECAUTIONS SHOULD**  
**BE STRICTLY FOLLOWED.**

**WARNING** - To reduce the risk of accidents:

- Only use this device for its intended purpose, as described in this manual.
- Before using Pearl, ensure that the skin to be treated is clean and dry and free from creams and oils.
- Do not use Pearl on the genitals, breasts or around the eyes.
- If you are pregnant, please consult your practitioner before using Pearl.
- Do not use Pearl on sunburned skin, open wounds, blisters or scars.
- If you have any skin condition other than those listed above, please consult your practitioner before using Pearl.
- If used correctly, Pearl should not cause any skin irritation. Nevertheless, it may cause irritation when used on sensitive skin. If this occurs, apply calming lotion after each treatment.
- Using Pearl on one area for too long may increase the possibility of skin irritation. If the unwanted hair has not been removed with three or four passes, wait for a few minutes before trying again on the same area.

**WARNINGS** - To reduce the risk of burns, fire, electric shock and injury:

- Do not place or store the device near a bath or sink.
- Do not place the device in or near water or other liquids.
- Do not attempt to touch a device that has fallen into water if it is plugged in and charging.
- This device must never be left unattended while plugged in.
- Do not use the head if the wire is twisted or damaged.
- Do not attempt to touch the wire inside the head. The plastic teeth of the head are designed to protect your skin during normal use but direct contact with the wire may cause injury.
- Switch off the device when changing the head.
- This device should only be used with the original charger supplied.
- Never use the charger if the cable or plug are damaged.
- Do not use Pearl near flammable objects.
- Always unplug the device straight after use and before cleaning or carrying out maintenance.
- Do not open the main unit. This may damage the device and void the warranty.
- Keep out of the reach of children.
- Do not use the device on pets.
- Do not attempt to repair Pearl if it becomes damaged. The product should be replaced by the manufacturer or distributor in your country.

Distributed by: **Windirect, SL**  
Rambla de Catalunya, 124,  
08008 Barcelona, SPAIN  
[www.windirect.es](http://www.windirect.es)

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