


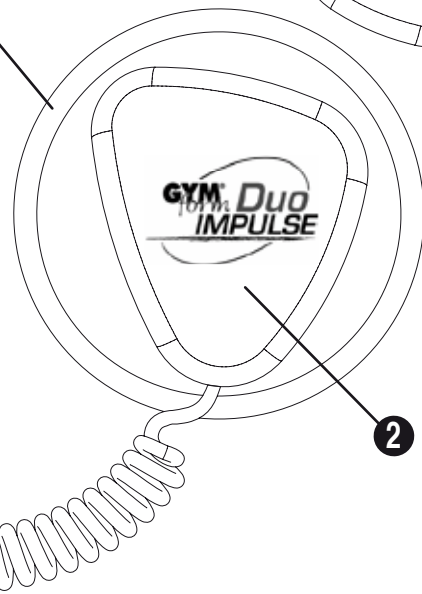
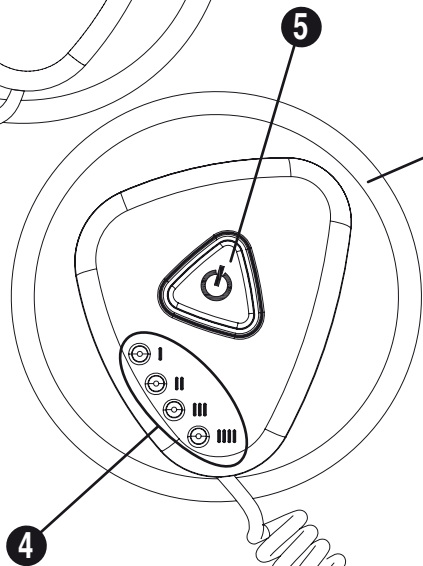
GYM[®] *form* Duo IMPULSE

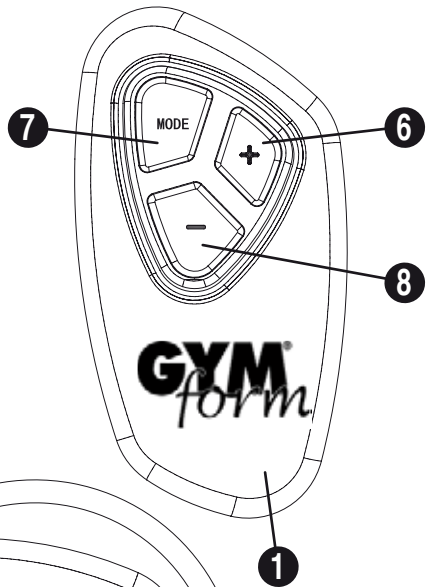- 🇬🇧 *Instruction Manual*
- 🇪🇸 *Manual de instrucciones*
- 🇫🇷 *Mode d'emploi*
- 🇩🇪 *Bedienungsanleitung*
- 🇮🇹 *Istruzioni per l'uso*
- 🇵🇹 *Manual de instruções*
- 🇳🇱 *Gebruiksaanwijzing*

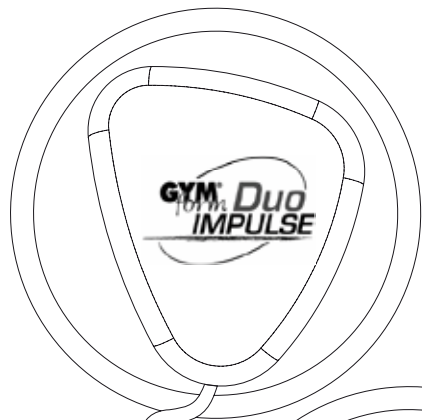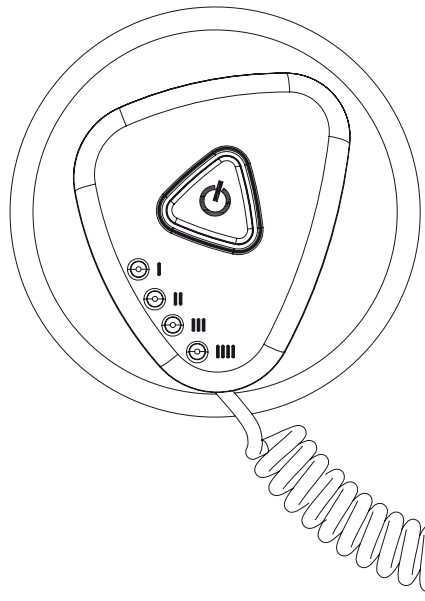GYMFORM® DUO IMPULSE INSTRUCTION MANUAL

Congratulations on your purchase of Gymform® DUO IMPULSE. Your safety is our primary concern, so please ensure that you familiarize yourself with all the information in this manual before using your Gymform® DUO IMPULSE. It is essential that anyone who intends to use this device reads and understands the instructions provided. An informed user is a safe user. Thank you for purchasing Gymform® DUO IMPULSE. We hope you enjoy your sessions.

WARNING

Before embarking on any exercise programme, including Gymform® DUO IMPULSE, consult your general practitioner and your physiotherapist.

Ensure that you read the instructions, warnings and precautions carefully before using your Gymform® DUO IMPULSE.

Gymform® DUO IMPULSE is designed for healthy adults. However, the device should be used with care. Read the following guidelines and the instruction manual before using the product.

- Applying electrodes in the thorax area may increase the risk of cardiac fibrillation.
- Using Gymform® DUO IMPULSE while connected to high-frequency surgical equipment may cause burns at the contact points of Gymform® DUO IMPULSE and may result in damage to the product.
- Operation of a short-wave therapeutic device or microwave oven in the general vicinity (e.g. 1 metre away) may cause instability in the output power of Gymform® DUO IMPULSE.
- This is a powerful electrical muscle stimulation device and may cause injury if you start with the most advanced exercise programme right away. Start with a lower level and gradually increase the intensity.
- Do not use on children or animals.
- Do not use while driving or sleeping or if you feel drowsy.
- Do not use in the bathroom or in other places with a high amount of moisture or humidity.
- This device may be used by children aged 8 and over, by people with reduced physical, sensory or mental abilities, and by people lacking in experience and knowledge, provided they are appropriately supervised or given suitable instructions on the safe use of the device and understand the dangers its use entails. Children should not be allowed to play with the device. The cleaning and maintenance to be performed by the user should not be done by children without supervision.
- Avoid using the device near equipment that emits electromagnetic radiation such as telecommunications antennas, microwave ovens and X-ray equipment, as it may interfere with their function.
- The life of the product depends on skin type, intensity of use and storage conditions. We recommend replacing them after approximately 30 uses to ensure maximum results. Replacement electrodes are available from the centre where you purchased the product. If you notice a premature decline in the electrodes' performance, moisten them with a wet rag.

DO NOT use Gymform® DUO IMPULSE:

- If you have an electronic implant (such as a pacemaker) or suffer from another heart problem.
- If you are, or think you might be, pregnant.
- If you have cancer, epilepsy or are under medical supervision for a cognitive disorder.
- If you are breastfeeding.
- If you have a high fever.

Women: before using your Gymform® DUO IMPULSE, you should wait:

- At least six months after the birth of your child (after consulting your doctor).
- A month after having an intrauterine device (IUD) fitted.
- At least three months after having a caesarean section (after consulting your doctor).
- Until after the first days of your menstrual period.

Before using Gymform® DUO IMPULSE, you should obtain permission from your doctor or physiotherapist:

- If you suffer from back problems. In this case, the device should only be used at low intensity levels.
- If you suffer from a serious disease or injury not mentioned in this section.
- If you were operated on recently.
- If you are diabetic and take insulin.

Do not place one receiver on your chest and another on your back in such a way that electricity passes through your heart. Do not place the receivers on genitals, scars, moles, varicose veins or the head.

- Do not place on recent scars, inflamed skin, infected areas or areas that are susceptible to acne, thrombosis or other vascular problems (e.g. varicose veins), or any part of your body with limited sensitivity.
- Do not place on areas with injuries or limited mobility (e.g. fractures or sprains).
- Do not place directly on metal implants. Place the receivers on the closest muscle.

Possible Adverse Reactions:

There have been isolated reports of skin reactions caused by the use of electro-stimulation devices, including allergies, prolonged skin redness and acne.

If any of these symptoms occur, switch off the device and consult your doctor.

If you have never exercised before or have not exercised for a long time: Consult your general practitioner.

Precautions:

- Do not use Gymform® DUO IMPULSE if your body is in contact with water.
- Do not use Gymform® DUO IMPULSE in combination with any other device that transmits electrical current to your body or with any other machine.
- Do not connect the metal snap fasteners or electrodes to any other object.
- Do not modify the product.
- Use only the electrodes supplied by Industex.

- Do not attach two self-adhesive electrodes together.
- If either of the electrodes is not attached properly, you will notice a tingling sensation.

SPECIAL FEATURES

- Patented dual wireless technology.
- Practical and compact microprocessor technology.
- Stylish, ergonomic design.
- Wireless units with oscillating dual-head technology.
- Electrodes with superadhesive conducting gel.
- 3 highly touch-sensitive control buttons (MODE, ON/- and OFF/-).
- 4 LED lights that indicate the selected programme.
- Each programme offers 10 intensity levels from beginner to advanced.
- 6 very advanced exercise programmes.
- Automatically switches off after 15 minutes so you do the recommended amount of exercise.
- Automatically switches off after 20 seconds if the device is turned on but not used. It also automatically switches off after 20 seconds if one of the electrodes is not in direct contact with the skin.

COMPONENTS:

- 2 Gymform® DUO IMPULSE receivers
- 4 self-adhesive electrodes
- 1 Gymform® DUO IMPULSE remote control
- 1 instruction manual

PRODUCT DESCRIPTION:

1. Remote control.
2. Gymform® DUO IMPULSE receivers.
3. Self-adhesive electrodes.
4. Receiver indicator lights.
5. Receiver ON/OFF button. 
6. ON/- button to turn on and increase intensity.
7. Mode button to change programme.
8. OFF/- button to turn off and reduce intensity.

BEFORE FIRST USE

1. To use Gymform® DUO IMPULSE, two 3V CR-2032 batteries are required for each receiver (2) and one is required for the remote control (1) (not included). Slide off the cover and insert the batteries with the positive end (+) facing up. Close the battery compartment cover.
2. Attach each receiver (2) to each self-adhesive electrode (3) using the two snap fasteners and make sure it is properly attached.
3. Remove the protective plastic backing from the self-adhesive electrode (3). (Keep the plastic backing to protect the electrodes after using them.)
4. Place the self-adhesive electrodes (3) on the parts of the body to be treated.

Parts of the Body to Work Out Using Gymform® DUO IMPULSE:

Each Gymform® DUO IMPULSE device allows you to work out 2 adjacent muscles at the same time:

1. Pectorals
2. Upper abs
3. Lower abs
4. Obliques
5. Glutes
6. Trapezius muscles
7. Calves
8. Inner thighs
9. Outer thighs

You can work out many different combinations of adjacent muscles. You choose the muscles you want to work out. Here are some suggestions:

10. Lower front thigh + upper front thigh
11. Upper outer front thigh + upper outer rear thigh
12. Arms: biceps + triceps
13. Arms: biceps + forearms

INSTRUCTIONS

After placing the self-adhesive electrodes (3) on the parts of the body to be treated, turn on the receivers (2) by pressing the button (5) with the  or ON/OFF symbol. Next, use the remote control (1) to choose from among 6 different programmes and 10 intensity levels.

- ON/+ button to turn on and increase intensity (6): Press this button to turn on the receivers and increase the stimulation intensity. Each time you press the button, the intensity level will be increased by one level (1-10). There are 10 levels in total.
- OFF/- button to turn off and reduce intensity (8): Press this button to turn off the receivers and reduce the stimulation intensity. Each time you press the button, the intensity level will be reduced by one level (1-10). There are 10 levels in total. The receivers can be turned off by pressing the button to the lowest intensity level.
- There are 10 stimulation intensity levels. Level 10 is the one that provides the most stimulation.
- You can change programmes at any time by pressing the Mode button (7). Each time you press it, a different programme will start at the minimum intensity level. You can change this to your desired intensity level using the + button.
- There are 6 Gymform® DUO IMPULSE programmes. Choose the one you like best.
- The maximum recommended time for use on each part of the body is 10 minutes.
- After using the product, switch off the receivers by pressing the ON/OFF button with the .

GYMFORM® DUO IMPULSE PROGRAMMES

Programme	Function	LED 1	LED 2	LED 3	LED 4
		Red 1	Yellow 1	Red 2	Yellow 2
1	Massage	On	Off	Off	Off
2	Muscle Building	Off	On	Off	Off
3	Anti-Cellulite Slimming I	On	On	Off	Off
4	Slimming II	Off	Off	On	Off
5	Anti-Cellulite	Off	Off	Off	On
6	Fitness	Off	Off	On	On

TECHNICAL SPECIFICATIONS

Batteries: the unit runs on 5x3V CR-2032

Maximum power consumption: 15 mA (RMS) (Load: 500Ω)

Impulse frequency: 0.2Hz ~ 250Hz

Automatic shutoff time: 15 ± 3 minutes

Modes: 6 programmes to choose from

Impulse intensity control: 10 intensity levels to choose from
Impulse intensity upon choosing a programme: the intensity level starts at one each time you select a different programme.

Dimensions: 24cm x 6cm x 2cm // 7cm x 3,5cm x 1,3cm

Weight: 140gr

Fig. 1

This symbol indicates that the device is protected against electrostatic discharge by BF-type double insulation.

ISO 7000-1641 User Guide

CLEANING AND MAINTENANCE

DO NOT WASH Gymform® DUO IMPULSE. WIPE WITH A DRY OR DAMP CLOTH.

Remove the batteries from your Gymform® DUO IMPULSE if you do not plan to use it for an extended period of time.

MAINTENANCE AND STERILIZATION OF ADHESIVE ELECTRODES

- Do not immerse in water.
- Keep in a cool, dry place at room temperature.
- Do not remove the insulating label.
- After using the device, always replace the protective plastic backings on the adhesive electrodes.
- Store in the bag supplied.
- Sterilize by wiping with wet fingers.

PRECAUTIONS WITH BATTERIES

Ensure the batteries are properly installed before turning the unit on.

Turn off the unit to change the batteries.

Remove the batteries if you are not going to use the device for an extended period of time.

Take the necessary precautions when disposing of the batteries. Do not expose them to fire.

DISPOSAL OF MATERIALS

A symbol of a crossed-out wheeled bin means you should find out about and follow local regulations about disposing of this kind of product.

Do not dispose of this product as you would other household waste. Dispose of this device in accordance with the corresponding local regulations.

Electrical and electronic devices contain hazardous substances that can have harmful effects on the environment and/or human health and should be recycled properly.

The battery should be disposed of separately. To protect the environment, all batteries, including rechargeable ones, should not be disposed of as you would other household waste. Take them to a disposal centre and ensure it complies with applicable legislation on battery disposal.

Quality Warranty:

This product is covered by a warranty against manufacturing defects subject to the time periods stipulated by the legislation in force in each country. This warranty does not cover damages resulting from inadequate use, negligent commercial use, abnormal wear and tear, accidents or improper handling.

Australia & New Zealand only

This product has the benefit of certain Consumer Guarantees. These are prescribed by the Australian Consumer Law & the New Zealand Consumer Guarantees Act 1993 both of which provide protection for consumers. There is no express warranty for this product in Australia or New Zealand. The above paragraph refers to other countries.

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1. Para usar Gymform DUO IMPULSE © se necesitan 2 pila CR-2032 de 3 V para cada receptor (2) y 1 para el mando a distancia (1) (no incluidas). Deslice la tapa e introduzca las pilas con la cara positiva (+) hacia arriba. Cierre la tapa del compartimento.
ISL GmbH
Best Direct International PTY LTD
Best Direct B.V.
ISL ITALY S.r.l.
Best Direct International Ltd,